


The 60 Second Denver Nasal Splint®

- Stabilization of the nose after surgery
- Reduction of edema after nasal surgery or trauma
- Moldable aluminum retains nasal bones in position
- Layers give cushioned pressure on subcutaneous tissue
- Dorsal pad reduces dead space and skin wrinkling
- Dorsal pad reduces possibility of polybeak
- Packaged non-sterile with sterile components



Splinting Steps

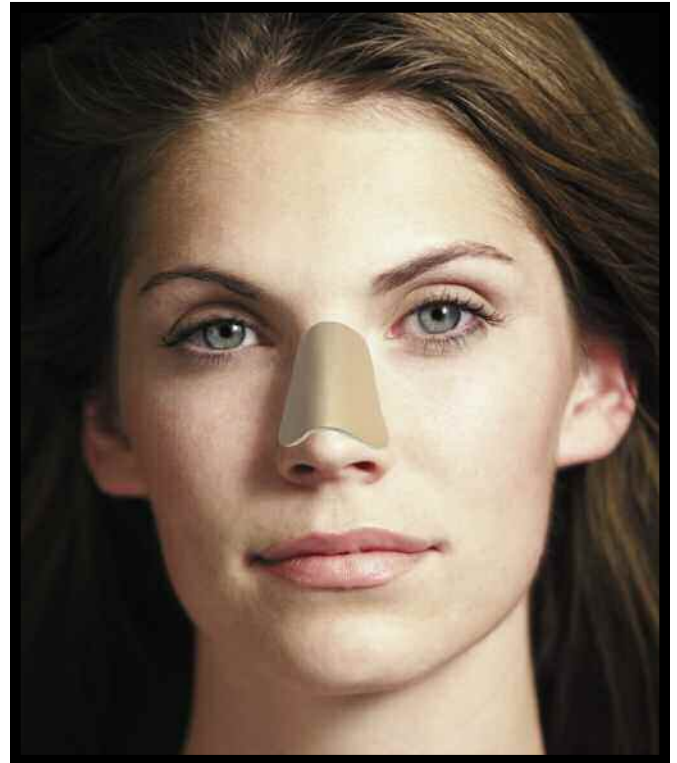

1- Alcohol and Skin Prep



2- Protecto Tape



3- Dorsal Pad



4- Adhesive Velcro




5- Aluminum Splint

- 60 seconds to apply
- Inconspicuous
- No operating room mess
- Water resistant
- One size per box



KIT COMPONENTS:

- Instruction sheet
- One alcohol wipe
- One skin prep
- One dorsal pad
- Pre-cut Protecto™ Strips
- One Denver adhesive Nasal Splint®

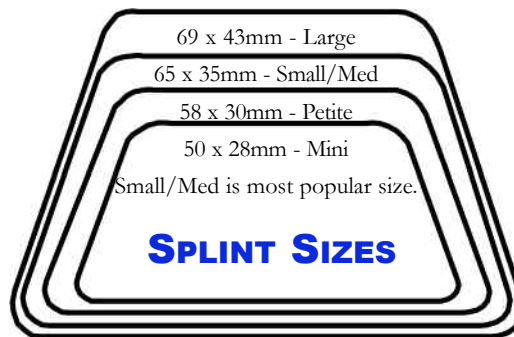

DESIGNER COLORS

Designer colors are available in all Denver Splint® models at no extra charge. One color per box. 10D-DS-series