



## RHINOCELL® SINUS PACKING

### Description

RHINOCELL® Sinus Packs are constructed of a unique sponge material using an exclusive formula. The small cell structure of RHINOCELL® together with a non-stick surface treatment provides a smooth surface, which maximizes comfort while minimizing tissue ingrowth. The high tensile strength of the sponge material prevents shedding or tearing. The packs are supplied compressed for ease of insertion and will expand as fluid is introduced, giving time for accurate positioning.

### Indications for use


For use as sinus packing.

### Supplied


Sterile

20 packs per dispenser carton


(2 packs per pouch)

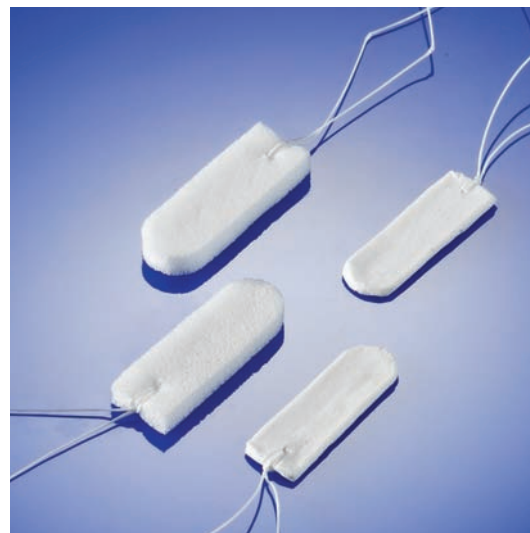

SP-350612

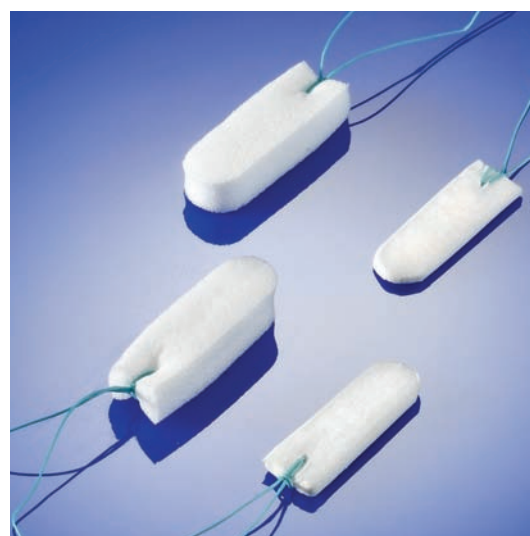

SP-350912

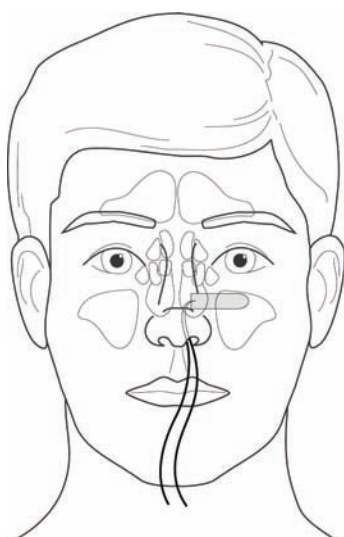

### Product Codes & Dimensions

REF	Rhinocell® Sinus Pack
SP-350612	3.5 x 0.6 x 1.2 cm
SP-350912	3.5 x 0.9 x 1.2 cm

### FEATURES

- Biocompatible lint and fiber-free sponge material
- Superior liquid absorption and wicking
- Supplied compressed for ease of insertion
- Extremely soft and pliable cell structure
- Attached lanyard to secure packing